



QAQC LAB

www.qclabequipment.com

Disposable Powder Sampler Large Volume, single use powder sampler


The disposable Powder Sampler is ideal for taking a single large volume sample from a single known point.

- Assembled & Packed in a Cleanroom
- Individually Bagged
- FDA & EU 10/2011 Conforming Materials
- BSE/TSE Free
- Full Batch Traceability
- ATEX Rated
- Available Sterile (gamma irradiated)



The disposable Powder Sampler is quick and easy to use:

1. Push the sampler into the product
2. Rotate handle to open up slots
3. After leaving the sampler for several seconds to allow powder to enter, rotate handle to close the sampler.
4. Remove sampler
5. Deposit the sample into a suitable container

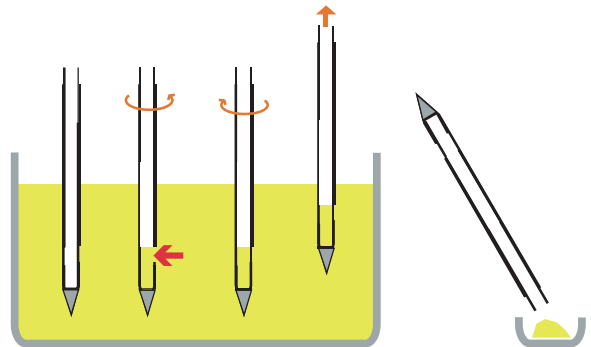

| Stock No. | Length | Vol. | Pre-Sterilized | Diameter | Material | Per Box |
|----------------|--------|-------|----------------|----------|----------|---------|
| 675 8240H1000 | 1000mm | 100ml | No | 40mm | HDPE | 20 |
| 675 8240H1000S | 1000mm | 100ml | Yes | 40mm | HDPE | 20 |

[QUOTE REQUEST](#) ✓

[CLICK FOR PRICING](#) ✓